




## FUME HOODS WITH EXTERNAL CONNECTION



▲ BH09D+ + BB09

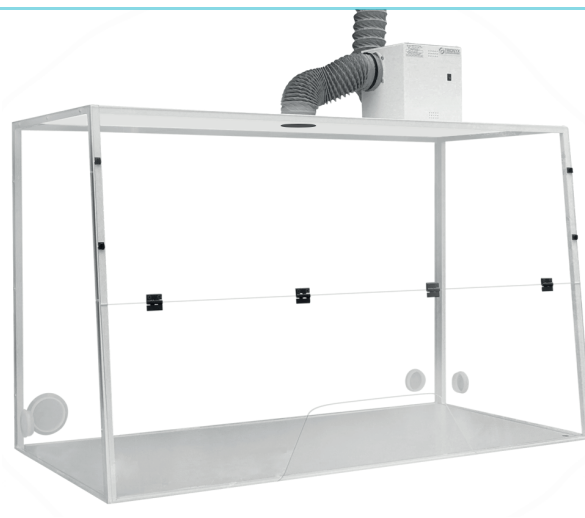

▲ BH15Z+ + BB15 + CDV-A + 2 KRC

### USE

These fume hoods are designed to be connected to a ventilation system(1) for the discharge of toxic and noxious fumes to the external of the building (use without active charcoal filter). It is imperative to add a ventilation box - model CDV-A (to be completed with 2 connection kits reference KRC or to be connected to a suitable extraction system).

### ADVANTAGES

- Hoods delivered ready to use (no assembly), to be connected only to the electrical network and to a fume extraction system.
- Transparent walls on all sides for optimal lighting of the work surface and immediate visualization of the products handled.
- Ergonomically designed working opening for safe and easy handling in the enclosure.
- Covers for the introduction of power and fluid supply cables.

### TECHNICAL CHARACTERISTICS

- 15/10th steel construction.
- White epoxy paint RAL 9010, highly resistant to chemicals.
- 4 sides made of 6 mm thick glass (PMMA)
- 4 covers for power supply passage.

### ADDITIONAL ACCESSORIES

- Ventilation box.
- Connection kit.
- Melamine worktop.
- Tempered glass worktop.
- Removable containment sump.
- External lighting.
- Fixed and rolling tables.

P/N	BH06D+	BH06Z+	BH09D+	BH09Z+	BH12D+	BH12Z+	BH15D+	BH15Z+	BH15G+	BH18Z+	BH18G+
Type of working opening	2 hands	Trapeze	2 hands	Trapeze	2 hands	Trapeze	2 hands	Trapeze	Large width	Trapeze	Large width
External dimensions W x D x H (mm)	600 x 650 x 745		900 x 750 x 945		1200 x 750 x 945		1500 x 750 x 945		1800 x 750 x 1240		
Internal dimensions W x D x H (mm)	545 x 505 x 735		845 x 700 x 935		1145 x 700 x 935		1445 x 700 x 935		1745 x 700 x 935		
Weight (kg)	23		34		44		54		70		
<b>Additional accessories</b>											
Containment sump + melamine worktop	BB06 (10 kg)		BB09 (11 kg)		BB12 (12 kg)		BB15 (13 kg)		BB18 (15 kg)		
Containment sump + laminated glass worktop	BV06 (15 kg)		BV09 (16 kg)		BV12 (17 kg)		BV15 (18 kg)		BV18 (20 kg)		
Ventilation box	CDV-A (8 kg)										
Connection kit	2 x KRC (2 kg)										
External side lighting	LUD (Exterior light with deflector - left hand side)										

(1) No model, even with an external connection, complies with the NF X 15-211 standard - (05/2009).