

EN
14470-2

EN
1363-1

30 APPROVED
EN14470-1

MADE IN
EUROPE



Fume hoods, filtering
cabinets - Ventilation

Lithium-ion
cabinets

**Flammable
cabinets**

Corrosive
cabinets

Toxic cabinets
and Multirisk

Containment
and cans

File cabinets
Anti-fire equipments

Showers and first-
aid equipments

30 MINUTES SAFETY CABINETS EN 14470-2 FOR GAS CYLINDERS



▲ 7635BG + RH5BG + EG1



▲ 7634BG + RHBG



Standardized labeling



1

Mounting rails with tension belt RH4BG / RH5BG (optional)




2

Shelf for one small gas cylinder P/N EG1 (optional)



3

Key locking system

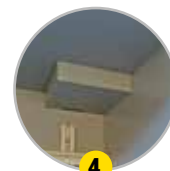

4

Vent with thermo-fuse

COMPLIANCE

- Tested and approved to EN 14470-2 (04/2006) and EN 1363-1 (06/2000).
- Standardized warning labels in accordance to EN 14470-2, ISO 3864, ISO 7010 and NF X08.003, and with Directive 92/58/EEC.

ADVANTAGES

- 30 minutes fire resistance according to the ISO 834 standard curve.
- Fire tested by a French test laboratory accredited to EN:17025 by COFRAC, according to the ISO 834 standard curve.
- Base with cover for easy movement by pallet truck, when empty.

PASSIVE SAFETY

- Double insulated panel construction offering high resistance to chemical attack.
- Light grey RAL 7035 frame with yellow RAL 1023 plated fields.
- Standardized symbol signage in accordance to EN 14470-2.

ACTIVE SAFETY

- Ventilation holes Ø 100 mm for possible connection.
- Thermo-expanding door seals.
- Ventilation ducts with hot-melt system to isolate the contents of the enclosure in the event of fire.
- Adjustable leveling feet
- Lockable doors.

OPTIONS

- To improve the longevity and safety of your cabinet, we recommend that you connect it to a forced ventilation system. For more information, please contact us.
- In order to increase the safety of your cabinet and to control the internal temperature, we propose a refrigeration system. Please contact us for more information.

All models are delivered without interior equipment. They must therefore be equipped with rails for the support of gas bottles, access ramps, etc.


RANGE 7630BG

EN
14470-2


EN
1363-1

30
APPROVED
EN14470-1


MADE IN
EUROPE

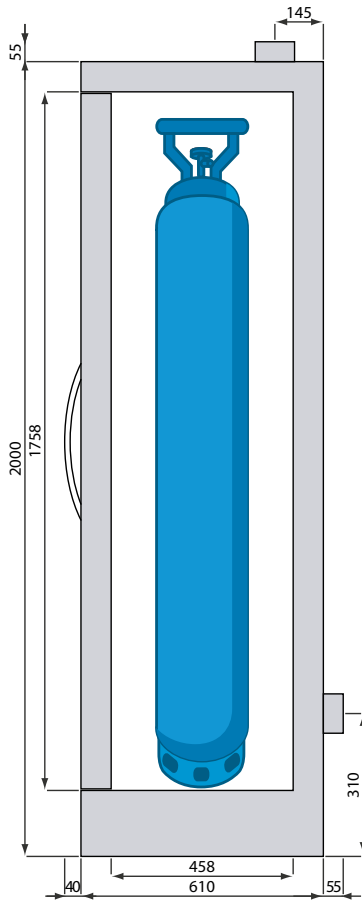

▲ 7634BG / 7635BG

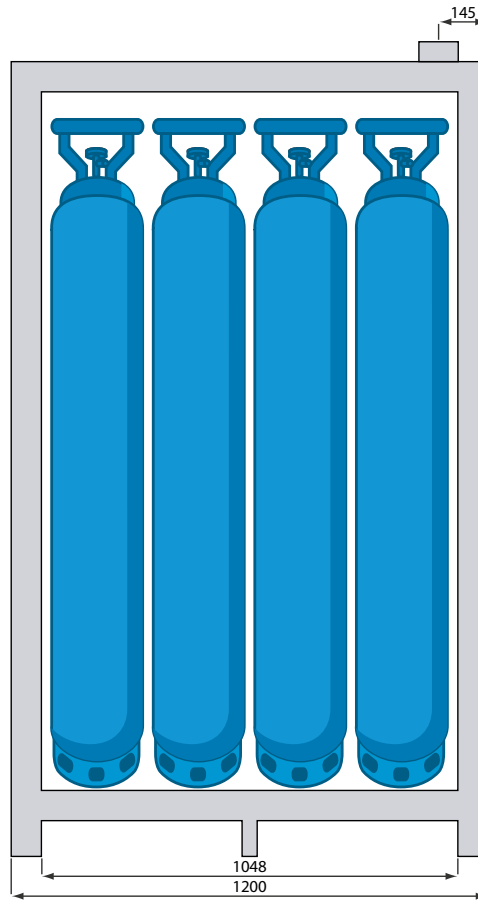

▲ 7635BG

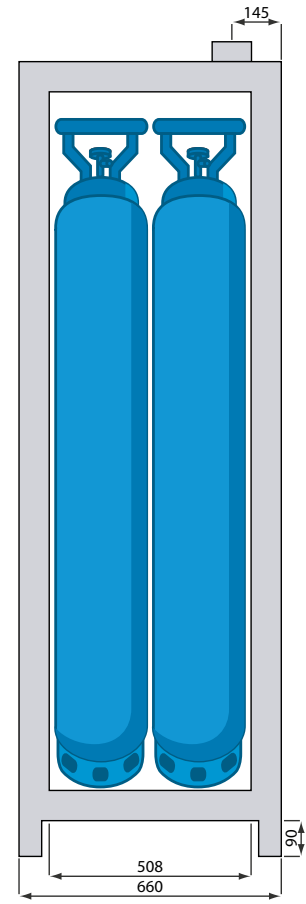

▲ 7634BG

P/N	Description	External dimensions H x W x D (mm)	Internal dimensions H x W x D (mm)	Bottle storage storage in B50 bottles	Recommended equipment	Weight (kg)
7634BG	Tall cabinet 1 door	2055 x 660 x 665	1753 x 508 x 458	2	1	246
RH4BG	Horizontal rails with straps for fixing 2 B50 bottles	-	-	-	1	2
7635BG	Tall cabinet 2 doors	2055 x 1200 x 665	1753 x 1048 x 458	4	1	447
RH5BG	Horizontal rails with straps for fixing 4 B50 bottles	-	-	-	1	3
RABG	Galvanized steel ramp for 7634BG or 7635BG	38 x 450 x 398	-	-	1	10
EG1	Additional shelf for 1 gas cylinder type B2, B5, B10, B11, B15 or B20	245 x 245 x 245	-	-	-	3

Fume hoods, filtering
cabinets - Ventilation

Lithium-ion
cabinets

Flammable
cabinets

Corrosive
cabinets

© ECOSAFE reserves the right to discontinue or change specifications or designs at any time without notice and without incurring obligations. 09/2021

Toxic cabinets
and Multirisk

Containment
and cans

File cabinets
Anti-fire equipments

Showers and first-
aid equipments